



BRILLIANT FUTURES™
MYOPIA MANAGEMENT PROGRAM



MiSight® 1 day is the innovative lens that goes beyond vision correction

MiSight® 1 day is the first and only contact lens FDA-approved to slow the progression of myopia* in children age 8-12 at initiation of treatment.**

CLINICALLY PROVEN OUTCOMES:

59% average reduction in SERE with MiSight® 1 day lenses^{1**}

52% average reduction in axial lengthening with MiSight® 1 day lenses^{1**}

9 out of 10 children wearing MiSight® 1 day lenses report seeing well while doing schoolwork, while playing outdoors, and while looking at a computer screen/playing video games.^{2††}



The MiSight® 1 day lens is the cornerstone of the Brilliant Futures™ Myopia Management Program, a comprehensive myopia management approach offered by CooperVision®.

***Indications for use:** MiSight® 1 day (omafilcon A) soft (hydrophilic) contact lenses for daily wear are indicated for the correction of myopic ametropia and for slowing the progression of myopia in children with non-diseased eyes, who at the initiation of treatment are 8-12 years of age and have a refraction of -0.75 to -4.00 diopters (spherical equivalent) with ≤ 0.75 diopters of astigmatism. The lens is to be discarded after each removal.

** Compared to a single vision 1 day lens over a 3 year period.

† From 1 month through 3 year visits


‡ From 1 week through 3 year visits

REFERENCES: **1.** Chamberlain P, et al. A 3-year Randomized Clinical Trial of MiSight Lenses for Myopia Control. *Optom Vis Sci.* 2019;96(8):556-567. **2.** Sulley A et al, Wearer experience and subjective responses with dual focus compared to spherical, single vision soft contact lenses in children during a 3-year clinical trial. AAO 2019 Poster Presentation.

HOW MISIGHT® 1 DAY WORKS:

The ActivControl™ technology in MiSight® 1 day uses vision correction zones and treatment zones within the lenses to slow axial elongation.¹⁻³

By including both types of zones in the lens, it simultaneously corrects the child's vision today, while creating myopic defocus in all gazes.¹

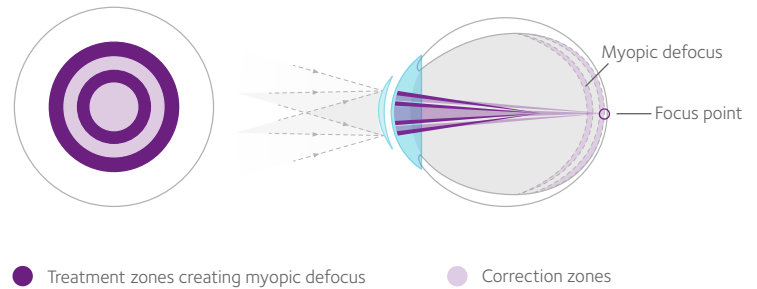

PARAMETERS:

Material	omafilcon A
Water content	60%
Base curve (mm)	8.7
Diameter (mm)	14.2
Sphere power	-0.5D to -6.00D (0.25D steps)
Optical design	ActivControl™ technology
Dk/t (at -3.00D)	28
Pack size	30
Replacement schedule	Daily disposable

Indications for use: MiSight® 1 day (omafilcon A) soft (hydrophilic) contact lenses for daily wear are indicated for the correction of myopic ametropia and for slowing the progression of myopia in children with non-diseased eyes, who at the initiation of treatment are 8–12 years of age and have a refraction of -0.75 to -4.00 diopters (spherical equivalent) with ≤ 0.75 diopters of astigmatism. The lens is to be discarded after each removal.

Talk to your CooperVision® Sales Representative to learn more about the Brilliant Futures™ Myopia Management Program, or visit [CooperVision.com/practitioner/myopia-management](https://www.CooperVision.com/practitioner/myopia-management)

